tover[®]


Tovertafel 3

June 2025

English	03
Français	09
Deutsch	15
日本語	21
Italiano	27
Nederlands	33
Svenska	39
Dansk	44
Norsk	51
Suomi	57
Español	63
Català	69
Galego	75
Euskara	81
Calibration marker	87

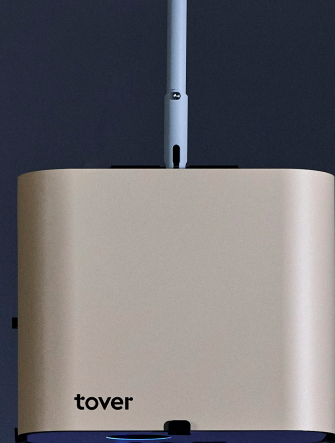Tovertafel 3 – User manual

Congratulations!

Well done on choosing Tovertafel. You now belong to a growing group of people who together bring fun and meaningful interaction to health care. The Tovertafel is not an ordinary product, but something special, developed together with health care professionals and target groups.

This manual provides a brief overview of the most important functionalities. For a detailed explanation, additional documentation and videos, please visit tover.care/ manual.

Do you want to get everything out of the Tovertafel? Then go to my.tover.care for tips, guidance and additional information. For example, here you can upload your own photos to the Tovertafel using ‘Lab’. This will make your games even more personal!

Purpose of the product

The Tovertafel stimulates movement, cognition and social interaction. The games are specially developed for people with dementia, an intellectual disability or other specific needs and are fun for everyone.

Because it's all about playing together!

Safety instructions

Read the manual and the safety instructions below carefully before use.

General warnings

- i** Safe use is only guaranteed if properly installed according to the manual.
- i** To hang the Tovertafel, use only items supplied by.
- i** The installer is responsible for damage caused by incorrect assembly.
- i** Only use the Tovertafel only as intended. Incorrect use invalidates the warranty.
- i** Do not hang anything from the appliance or hanging system.
- i** Do not lean or push against the device.
- i** Only use the appliance indoors, in dry and dust-free areas.
- i** Never open the appliance or insert metal objects.
- i** Do not touch or look into the lens for long periods of time.
- i** Keep children away during installation.
- i** The Tovertafel is not a medical device.
- i** Snap the handle firmly into the hook and check that it locks.
- i** Do not drop the appliance.

Electricity

- i** Only connect to an electrical outlet that is earthed and switched off.
- i** Turn off the appliance switch on the Tovertafel before you install or remove the plug.
- i** Make sure the cord is long enough.
- i** Do not expose to moisture, naked flames, heat, or sunlight.
- i** Remote control batteries
- i** Only use the correct type of battery (CR2032).

- i** Observe the correct + and – marks when positioning.
- i** Remove the batteries when you do not use the remote control for a long time.
- i** Do not expose to moisture, fire, heat, or sunlight.
- i** Keep out of reach of children.
- i** Immediately call for medical assistance in case of ingestion or when coming into contact with leaking fluid.
- i** Dispose of batteries according to local regulations.

Contents of the box

- Tovertafel 3
- Remote control + battery
- User manual
- Calibration marker

Delivered separately:

- Ceiling bracket + appliance plug
- Spare remote control
- Playing tiles (exclusive to the “Including Education” license)

Technical specifications

- Resolution: 1280 × 800 px
- Light source: 4 × LED
- Dimensions: 41 × 23 × 36 cm
- Weight: ±10.5 kg
- Power supply: 100–240 Vac / 50–60 Hz
- Power: 190 W
- Connectivity: Wi-Fi, ethernet, 4G, USB


Installation

Assembly

The Tovertafel can be installed on the ceiling bracket or a mobile stand (sold separately). Both products have their own manual. With multiple ceiling brackets, you can hang the Tovertafel in different locations. Additional ceiling brackets are available separately.

Parts of the Tovertafel

1. Power socket
2. Device switch
3. Handle
4. Ethernet port
5. USB ports + indicator LED
6. Infrared light source
7. Sensors
8. Projector lens

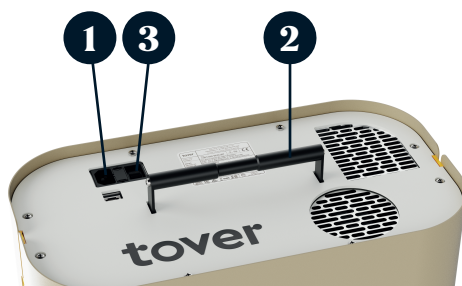Turn on the Tovertafel

1. Plug the power cord into an earthed wall socket.
2. Turn on the device switch. The indicator LED turns on and flashes red.
3. Wait until the indicator LED lights up red continuously.
4. Press the power button on the remote control.
5. The indicator LED flashes green.
6. After a while, the indicator LED lights up green continuously. The projection appears and the system is started.

The Tovertafel is now ready for use.

Turn off the Tovertafel


1. Press the power button on the remote control.
2. The indicator LED flashes red.
3. Wait until the LED is solid red.

The system is now safely shut down. If necessary, it is now possible to turn off the appliance switch and safely disconnect the plug.

Get started

Remote control

1. On/off
2. Back
3. Start/pause
4. Volume
5. Sound on/off
6. Navigation menu
7. OK
8. Menu
9. Info

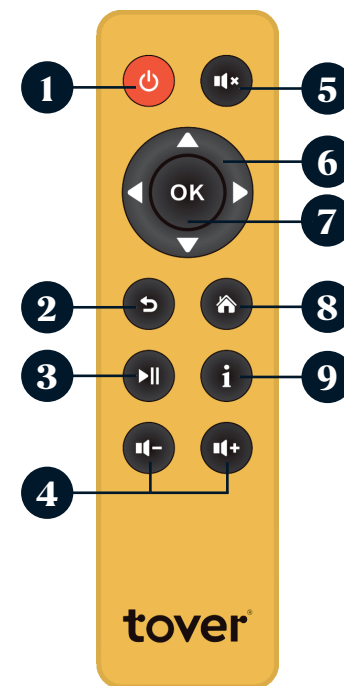Menu overview

Access the menu using the ‘Menu’ button on the remote control.

Depending on the activated license, the Tovertafel menu may change slightly. Use the arrow keys on the remote control to navigate through the menu. At the very top you will find the following tabs:

1. All games
2. Collections (optional)
3. Use (optional)
4. Settings

All games

If you have multiple licenses, the name of the active licenses will appear below the main menu. You can switch packages by pressing right or left.

Below are the available games. The games are categorized at the game level. Each row starts with a 'shuffle' symbol. Selecting this option will automatically randomly play the games from that row.

Collections (optional)

Here you will find specially created collections that make it easier to choose the right experience. Collections are composed based on theme, purpose or functionality (such as high-contrast games).

Insights (optional)

This menu provides insight into the use of the Tovertafel.

Settings > Calibrate playing surface

Use this menu to resize the projection to the table. Follow the on-screen instructions.

Settings > Set Table Color

Set the color of your table here so that the sensor works optimally. The options are 'dark' and 'light'.

Settings > Calibrate Sensor

If the games do not respond properly, you can reset the sensor here. Use the 'Auto Calibration'. If this does not work properly, you can calibrate manually using the paper calibration marker.

Settings > Network

The Tovertafel 3 comes standard with a mobile data connection for service, activation and critical updates. Major software updates are only offered over Wi-Fi or Ethernet connection. We recommend connecting the Tovertafel with an Ethernet cable, or a Wi-Fi point. To do this, go to 'Settings > Network' and follow the instructions in the menu.

Caution:

- *There are no additional charges for the mobile data connection and if desired it is possible to turn off the connection via the network menu.*
- *Open WiFi networks, networks with a login page, or with strict security settings are not usually suitable for use with the Tovertafel.*

Set language

From this menu you can choose from an extensive list of languages.

Info

In this menu you will find information about your Tovertafel, such as the serial number and the software version.

Tips and inspiration

Visit my.tover.care for inspiration, how-to videos, valuable tips and detailed information. Here you can also find 'Tover Lab'. With 'Tover Lab' it is possible to upload your own images to your Tovertafel.

Software and Licenses

Updates are automatically installed when connected to the Internet. It is possible to update the Tovertafel offline, using a USB stick. Check out my.to.care/service-center/updates for this. The update menu will tell you if an update is available.

Licenses can be purchased or renewed through your Tovertafel distributor. A renewal reminder will automatically appear in the menu.

Purchased licenses can also be activated offline using a code in the 'Unlock subscription' menu. This menu appears automatically when a Tovertafel is offline. Open the menu and follow the instructions.

Other

Maintenance

- Clean the Tovertafel and Remote Control regularly with a dry cloth.
- Over time, dust can get on and into the housing. You can request a service from your distributor.

Waste disposal

- The Tovertafel and Remote Control are electronic devices. Do not dispose of them in household waste.
- Remove the battery from the remote control and dispose of it separately.
- Follow local guidelines for handling electronic waste and batteries.


Regulations

The Tovertafel 3 complies with the applicable regulations in the EU, USA, Canada, Australia/New Zealand and Japan.

For full certifications and legal statements: tover.care/compliance

Support

For contact details, manuals, videos and help:

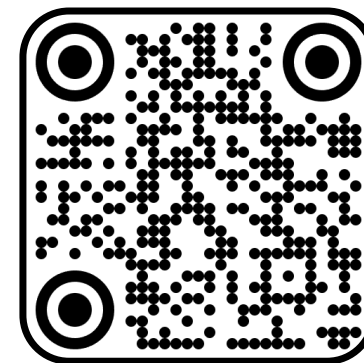Warranty conditions Tovertafel 3

Congratulations on purchasing the Tovertafel 3! We are committed to ensuring your satisfaction with the product and providing you with the best possible service. Please keep this document for your records.

Warranty

The Tovertafel 3 and its accessories are guaranteed to be free from defects in materials and workmanship under normal use and maintenance. This includes periodic software updates and support for licensed games.

Normal use means: Use according to the instructions in the manual, including recommended maintenance and software updates.

Period

The standard warranty period is 24 months from the date of activation of the software license.

The warranty period will not be extended after a repair or replacement.

New or functionally equivalent parts may be used for repairs or replacements.

Terms and Conditions

This warranty applies only to the first purchaser and is not transferable. The Tovertafel must be purchased for direct use within the organization, and not for resale.

Submit Warranty Claim

Please contact your local distributor, who will assist you in resolving the issue and submitting the warranty claim.

An overview of distributors can be found at tover.care/manual.

If necessary, the distributor will contact the service team.

Replacement device: In some markets, a temporary replacement table is available, which may be subject to a fee. Please check with your distributor for the terms and conditions.

Excluded

The following situations are excluded from warranty:

- Damage due to misuse, accident, alteration, incorrect installation or external influences (such as overvoltage, fire or water damage).
- Cosmetic wear, packaging damage, or wear of parts such as the projector lamp (up to 20,000 hours after normal use).
- Software other than the licensed games of Tover.
- Use outside of the recommended environmental conditions or with accessories that are not approved.

Regulations

This warranty replaces, to the extent permitted by law, all other legal warranties. This warranty does not affect legal consumer rights that may apply by country or region.

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Tover Health Tech B.V.